


μ-Slide 18 Well - Flat

Instruction Manual

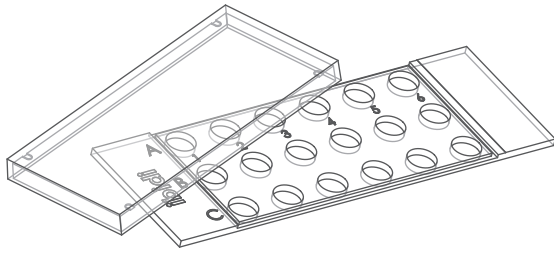

The ibidi labware is comprised of a variety of μ-Slides, μ-Dishes, and μ-Plates, which have all been designed for high-end microscopic analysis of fixed or living cells. The high optical quality of the ibidi Polymer Coverslip is similar to that of glass, enabling a variety of microscopy techniques with uncompromised resolution and choice of wavelength.

The μ-Slide 18 Well - Flat is a convenient multiwell plate for quick and easy testing of experimental parameters such as antibody dilutions, cell seeding densities, or optimal drug concentrations. Unlike standard 96-well plates, it enables full well-area observation using phase contrast microscopy, as the flat design minimizes disturbing refraction of the liquid surface meniscus.

This document applies to the following products:

81826	μ-Slide 18 Well - Flat ibiTreat
81821	μ-Slide 18 Well - Flat Uncoated

Material

The μ-Slide 18 Well - Flat is made of a polymer that has the highest optical quality. The ibidi Polymer Coverslip bottom exhibits extremely low birefringence and autofluorescence, similar to that of glass. It is not possible to detach the bottom from the upper part. The slide is intended for one-time use and is not autoclavable, since it is only temperature-stable up to 80°C/175°F. Please note that gas exchange between the medium and the incubator's atmosphere occurs partially through the polymer coverslip, which should not be covered.

Optical Properties of Polymer Coverslip

Refractive index (589 nm)	1.52
Abbe number	56
Thickness	No. 1.5 (180 μm)
Material	Polymer



WARNING – The ibidi Polymer Coverslip is compatible with certain types of immersion oil only. A list of suitable oils can be found in the Section "Immersion Oil".


Geometry

The μ-Slide 18 Well - Flat provides standard slide format according to ISO 8037/1.

The well-to-well distance of 9 mm (as in 96 well plates) allows the use of multichannel pipettes.

Specifications

Outer dimensions (w × l)	25.5 × 75.5 mm ²
Number of wells	18
Well diameter	5 mm
Volume per well	30 μl
Height with/without lid	5.0/1.6 mm
Growth area per well	0.2 cm ²
Coating area per well	0.25 cm ²
Bottom	ibidi Polymer Coverslip

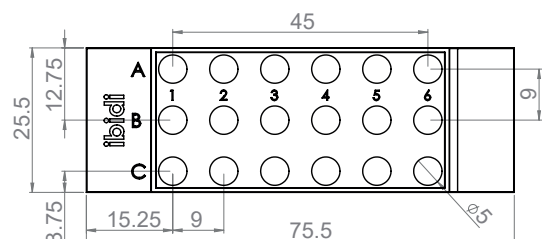

Shipping and Storage

This product is sterilized and sealed in a gas-permeable packaging. The shelf life under proper storage conditions (in a dry place, no direct sunlight) is outlined in the following table.

Conditions	
Shipping conditions	Ambient
Storage conditions	RT (15–25°C)

Shelf Life	
ibiTreat, Uncoated	36 months

Surface

The μ-Slide 18 Well - Flat is available with either an ibiTreat or an Uncoated surface.

The tissue culture-treated, hydrophilic ibiTreat surface of the ibidi Polymer Coverslip is ideal for culturing adherent cells. It ensures excellent cell adhesion without the necessity for any additional coatings. Nonetheless, extracellular matrix (ECM) protein coatings can be applied to the ibiTreat surface without any restrictions, if required.

The hydrophobic Uncoated surface of the ibidi Polymer Coverslip offers weak cell adhesion unless pre-coated with an ECM protein. You can apply coatings to the Uncoated surface without any restrictions. This surface is suitable for culturing adherent cells that require a specific coating.

For establishing a particular coating, we advise testing your procedure on both ibiTreat and Uncoated surfaces, as proteins and biomolecules may adhere differently to hydrophilic or hydrophobic surfaces.

Coating

Detailed information about coatings is provided in [Application Note 08: Coating Protocols for ibidi Labware](#).

In short, specific coatings are possible following this protocol:

1. Prepare your coating solution according to the manufacturer's specifications. Adjust the concentration to a coating area of 0.25 cm² and a volume of 30 μl per well.
2. Apply 30 μl per well and leave it at room temperature for at least 30 minutes.
3. Aspirate the solution and wash with the recommended protein dilution buffer.
4. The coated slide is ready to be used. Be aware that allowing the coated surface to dry out is not recommended, as some coating proteins may degrade upon drying.

Seeding Cells

1. Trypsinize and count the cells as usual. Dilute the cell suspension to the desired concentration. Depending on your cell type, application of a $0.7\text{--}1.7 \times 10^5$ cells/ml suspension should result in a confluent layer within 2–3 days.
2. Apply 30 μl cell suspension per well. Avoid shaking, as this will result in inhomogeneous cell distribution.
3. Cover the slide with the supplied lid. Incubate as usual (e.g., at 37°C and 5% CO₂).



TIP – The volume of a single well is very small. Depending on your cell type, the medium might be consumed after some hours. For incubations lasting longer than a few hours, we recommend refreshing the cell culture medium daily by aspirating and refilling.

Microscopy

To image your cells, no special preparations are necessary. Living or fixed cells can be directly observed, preferably on an inverted microscope. The bottom cannot be removed. For optimal results in fluorescence microscopy and for storage of fixed and stained samples, ibidi provides mounting media that are optimized for ibidi labware:

Cat. No. 50001: [ibidi Mounting Medium](#)

Cat. No. 50011: [ibidi Mounting Medium with DAPI](#)



TIP – When each well is properly filled with 30 µl of medium, the liquid surface is planar and in alignment with the μ-Slide surface. This enables clear, uniform phase contrast imaging across the entire well area.

Chemical Compatibility

The following table provides some basic information on the chemical and solvent compatibility of the μ-Slide 18 Well - Flat. For a full list of compatible solvents and more information on chemical compatibility, visit [ibidi.com/chemicals](https://www.ibidi.com/chemicals).

Chemical / Solvent	Compatibility
Methanol	Yes
Ethanol	Yes
Formaldehyde	Yes
Acetone	Yes, without lid
Mineral oil	No
Silicone oil	Yes
Immersion oil	See Section “Immersion Oil”

Immersion Oil



WARNING – When using oil immersion objectives with the ibidi Polymer Coverslip, use only the immersion oils specified in the table below. The use of any non-recommended oil could damage the ibidi Polymer Coverslip. The resulting leakage may harm objectives and microscope components. All immersion oils that are not listed in the table below should be considered as non-compatible.

Company	Product	Ordering No.	Lot Number	Test Date
ibidi	ibidi Immersion Oil 2	50102	24-07-04	07/2024
Cargille	Type A	16482	100592	01/2017
Cargille	Type HF	16245	92192	01/2017
Carl Roth	Immersion oil	X899.1	414220338	01/2017
Leica	Immersion Liquid	11513859	n.a.	03/2023
Leica	Immersion Liquid Type G	11513910	n.a.	04/2024
Nikon	Immersion Oil F2 30cc	MXA22192	n.a.	01/2020
Nikon	Silicone Immersion Oil 30cc	MXA22179	20191101	01/2020
Olympus	Silicone Immersion Oil	SIL300CS-30CC	N4190800	01/2017
Zeiss	Immersol 518 F	444960-0000	220211	03/2023
Zeiss	Immersol 518 F (30 °C)	444970-9010	220816	03/2023
Zeiss	Immersol 518 F (37 °C)	444970-9000	220302	03/2023
Zeiss	Immersol W 2010	444969-0000	101122	04/2012
Zeiss	Immersol Sil 406	444971-9000	80730	03/2023
Zeiss	Immersol G	462959-9901	211117	03/2023

For research use only!

Further information can be found at ibidi.com. For questions and suggestions, please contact us by e-mail at info@ibidi.com or by telephone at +49 (0)89/520 4617 0.
© ibidi GmbH, Lochhamer Schlag 11, 82166 Gräfelfing, Germany.